

## Red Light Therapy



### D: Original McLaren Red LED Torch - Usage Instructions

Photonic Therapy can provide a simple and safe healing process. To turn on the Original McLaren Red LED Torch simply turn the black grip to left or right.

Apply the Red Light to unbroken skin and ideally upon identified acupuncture points or over areas experiencing pain or stiffness. Please refer to Photonic Therapy Treatment Packs for more details (Human, Horse & Canine Kits).

### E: Product Support

Should you encounter problems with this device please follow these steps:

- Read the Manual
- Contact the WPT Support centre.

In Australia, support operates during normal business hours Monday – Friday.

- Telephone (0414) 515 400 [ask to speak with Trevor]
- Email [trevorandlaura@bigpond.com](mailto:trevorandlaura@bigpond.com)
- Website [www.wozenphotonicttherapy.com](http://www.wozenphotonicttherapy.com)

Trevor & Laura Wozencroft



## Red Light Therapy




### User Manual

## The Original McLaren Red LED Torch



*The Original McLaren Torch*

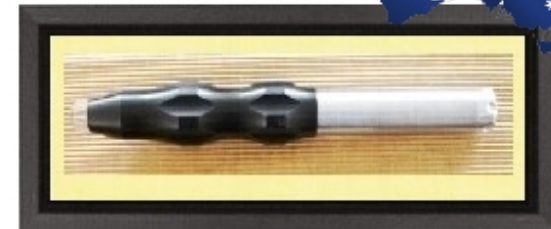

International Institute for  
Complementary Therapists  
**IICT MEMBER**  
[www.iict.com.au](http://www.iict.com.au)

# Red Light Therapy



## CONTENTS

	Page
A: Introduction	3
1.1 Package Contents	
1.2 Features & Specifications	3
B: LED Safety Information	3
C: Battery Installation & Replacement	3
D: Usage Instructions	
1.1 Dog / Cat / Humans Treatments	4
1.2 Medium / Large Animals (Humans / Horses)	4
E: Product Support	4



Made in Australia

Developed in Australia

## A: Introduction

Thank you for purchasing the **Original McLaren Red LED Torch**. WPT wishes you every success with your new LED treatment Light.

### 1.1 Package Contents

Carefully unpack and remove the contents. If any of the following items are missing or damaged, contact Wozenphotonicttherapy (WPT).

- **Original McLaren Red LED Torch**
- 2 x AA Batteries
- Optional Photonic Therapy Kit(s) [Human, Equine & Canine]
- Optional Sturdy Storage Case
- User Manual [Equine / Canine]
- Simple Health Maintenance Information [Humans]

### 1.2 Features & Specifications

#### The **Original McLaren Red LED Torch**

The 5 x Light Emitting Diodes (LED) output a 660 nm Hyper Red Light. The LEDs output up to  $3.9 \times 10^{17}$  photons / second / cm at 400 mW. Five year (5 YR) warranty for workmanship and product.

## B: LED Safety Information

The **Original McLaren Red LED Torch** is a **Class 1 / Low Risk LED** (battery – operated) **Torchlight**.

**WARNING** – Do not look directly at or shine Red LED light beam at eyes as this may cause brief visual impairment.

## C: Battery Installation & Replacement

To install / replace the two AA batteries provided simply unscrew Torchlight base to insert and replace to seal the Batteries. REMOVE batteries from LED Light if not used for prolonged time period to avoid corrosion.